ACCU-CHEK[®] *Performa*

BLOOD GLUCOSE METER / GLYCOMÈTRE

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Contents

Chapter 1: Introduction	7
The Accu-Chek Performa Meter	7
Coding the Meter	9
Adjusting the Time and Date – First Time Use.....	11
Using the Accu-Chek Performa System.....	13
Chapter 2: Patient Testing for Healthcare Professionals	15
Sample Handling	15
Performing a Blood Glucose Test	16
Flagging Test Results.....	19
Unusual Test Results.....	20
Recommending Alternate Site Testing to Patients	22
Chapter 3: Control Testing for Healthcare Professionals	23
Introduction.....	23
Testing Intervals.....	23
About the Control Solutions	24
Performing a Control Test	25
Control Test Results	29

Chapter 4: Linearity Testing for Healthcare Professionals	31
Introduction	31
Testing Intervals.....	31
Performing Linearity Tests	32
Linearity Test Results.....	33
Chapter 5: Meter Memory, Setup, and Downloading	35
Memory	35
Viewing Results	36
Meter Setup	37
Setting the Time and Date.....	38
Setting the Beeper On/Off	41
Setting the Alarm Clock Function	44
Setting the Hypoglycemic (Hypo) Alarm Function.....	48
Transferring Data Directly to a Computer or PDA Using Specialized Software and Infrared Cable	51
Chapter 6: Maintenance and Troubleshooting	53
Changing the Battery.....	53
Cleaning the Meter.....	54
Maintenance and Troubleshooting.....	55
Screen Messages and Troubleshooting	56

Chapter 7: Technical Information	63
Product Limitations	63
Specifications.....	63
Product Safety Information.....	66
Disposing of the Meter	66
Explanation of Symbols	67
Guarantee	68
Additional Supplies.....	68
Chapter 8: Blood Glucose Testing for Consumers	69
Performing a Blood Glucose Test	70
Flagging Test Results.....	74
Alternate Site Testing (AST).....	75
Unusual Test Results.....	80
Symptoms of High or Low Blood Glucose.....	82
Chapter 9: Control Testing for Consumers	83
Why Perform Control Tests.....	83
About the Control Solutions.....	84
Performing a Control Test	85
Understanding Control Test Results.....	89
Index	93

The Accu-Chek Performa System

The Accu-Chek Performa Meter is for quantitative blood glucose testing using Accu-Chek Performa Test Strips.

Intended Use

Test strips for quantitative blood glucose measurements in fresh capillary or arterial blood, as well as in anticoagulated venous blood. Refer to the test strip package insert for acceptable anticoagulants.

For use only with the Accu-Chek Performa Meter.

The meter may be used by people with diabetes to test their own blood glucose as well as for blood glucose determination by healthcare practitioners in a professional setting.

Suitable for self-testing

The system includes:

- **Accu-Chek Performa Meter with battery**
- **Accu-Chek Performa Test Strips and code chip**
- **Accu-Chek Performa Control Solutions**

Any object coming into contact with human blood is a potential source of infection (see: Clinical and Laboratory Standards Institute: Protection of Laboratory Workers from Occupationally Acquired Infections; Approved Guideline – Third Edition; CLSI document M29-A3, 2005).

Important Information About the Meter

- This meter is designed for testing fresh whole blood samples. The meter is for in vitro (outside the body) use. It should not be used to diagnose diabetes.
- The meter requires Accu-Chek Performa Test Strips. Other test strips will give inaccurate results.
- The meter comes with a preset time and date which may need to be set to the local time zone.

Patients receiving peritoneal dialysis using solutions containing icodextrin (e.g. Extraneal) should not use the Accu-Chek Performa Test Strip. It could cause an overestimation of glucose results.

Chapter 1: Introduction

The Accu-Chek Performa Meter

Display –

Shows results, messages, and results stored in memory.

Right and Left Arrow Buttons –

Press to enter memory, adjust settings, and scroll through results.

Test Strip Slot –

Insert test strip here.

Infrared (IR) Window –

Transfers data from the meter to a computer or PDA.

On/Off/Set Button –

Turns the meter on or off and sets options.

Battery Door –

Flip open the battery door by pushing the tab in the direction of the arrow.

Code Chip Slot –

Insert code chip into this opening.

Contacts –

Insert this end of the test strip into the meter.

Yellow Window –

Touch blood drop or control solution here.

**Test Strip****Test Strip Container****Control Solution Bottle**
(example)**On/Off/Set Button****Top View****Infrared (IR) Window****Code Chip**
(example)**Battery** – Insert with (+) symbol facing up.

Coding the Meter

1. Make sure the meter is off.

2. Turn the meter over.

3. Remove the old code chip (if there is one in the meter) and discard it.

Change the code chip every time a new box of test strips is opened!

4. Turn the code chip over so the code number faces down. Push the code chip into the code chip slot until it stops.



5. Leave the code chip in the meter until a new box of test strips is opened.

Notes:

- Do not force the code chip into the meter. It is designed to go into the meter only one way.
- If code “- - -” is displayed, insert a code chip into the meter.

Adjusting the Time and Date – First Time Use

The meter is preset with the time and date which may need to be adjusted to the local time zone. Setting the correct time and date in the meter is important to review the results in memory or to download the results to a computer.

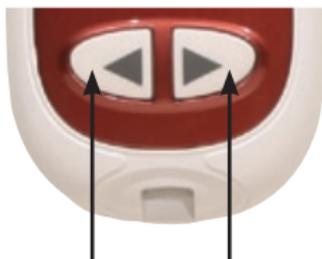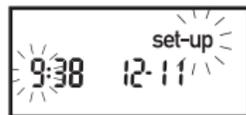decrease




increase

1. Press and release  to turn the meter on. The time and date appear on the display. "Set-up" and the hour flash. If the correct time and date are

2. displayed, press and **hold**  to exit. Press and release  or  to decrease or increase the hour. Press and **hold**  or 

3. down to scroll faster. Press and release  to set the hour. The minutes flash.

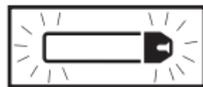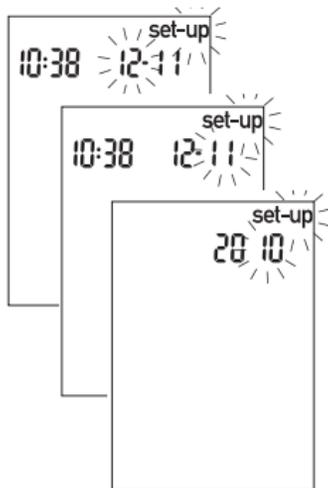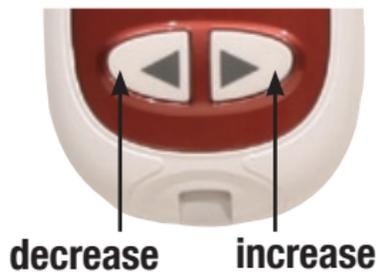4. Press and release  or  to adjust the minutes. Press and release  to set the minutes.

5. Repeat to set the day, month, and year.
6. Once the year is set, press and **hold**  until the flashing test strip symbol

appears. The time and date are now set.

7. See Chapter 5 “Meter Setup” to set other options.

Note: After the battery is changed, the meter automatically prompts the operator to confirm the time and date when it is turned on.

Using the Accu-Chek Performa System

- Only use Accu-Chek Performa Test Strips.
- Replace the code chip every time a new box of test strips is opened.
- Store the unused test strips in their original container with the cap closed.
- Close the container tightly immediately after removing a test strip to protect the test strips from humidity.
- Use the test strip immediately after removing it from the container.
- Check the expiry date on the test strip container. Do not use the test strips after that date.
- Store the test strip container and meter in a cool dry place.
- Refer to the test strip package insert for test strip storage and system operating conditions.
- Do not apply blood or control solution to the test strip before inserting it into the meter.

Do not store test strips in high heat and humidity! Heat and moisture can damage the test strips.

Consumers:

Please proceed to Chapters 5 through 9. Chapters 2, 3, and 4 are for healthcare professionals only.

Chapter 2: Patient Testing for Healthcare Professionals

Sample Handling

Always wear gloves and follow your facility's infection control procedures when handling blood-contaminated items. Always adhere to the recognized procedures for handling objects that are potentially contaminated with human material. Follow the hygiene and safety policy of your laboratory or institution. Prepare the selected blood collection site per facility policy.

A drop of fresh, whole blood is required to perform a blood glucose test. Fresh venous, capillary, or arterial blood may be used. Caution should be taken to clear arterial lines before blood is drawn. Blood glucose determination with venous or arterial blood must be performed within 30 minutes of sample collection. Avoid air bubbles with the use of pipettes. Refer to the test strip package insert for appropriate anticoagulants for fresh venous blood specimens.


Performing a Blood Glucose Test

16


Before performing a blood glucose test on patients, set up the meter correctly. The meter, a test strip, and a disposable lancet or blood collection device are required.

1. Prepare the lancet or blood collection device.
2. Insert the test strip into the meter in the direction of the arrows. The meter turns on.
3. Make sure the code number on the display matches the code number on the test strip container. If the code number is overlooked, remove the test strip and reinsert it into the meter.

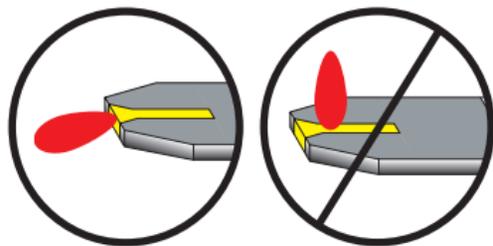4. A test strip and flashing blood drop symbol appear on the display.
5. Obtain blood from the patient per facility policy.
6. Touch the blood drop to the **front edge** of the yellow window of the test strip. Do not put blood on top of the test strip. When  flashes, sufficient blood is in the test strip.

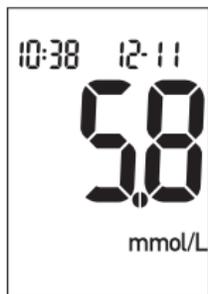
7. The result appears on the display. To flag a result, leave the test strip in the meter. (See the next section.) Otherwise, remove and discard the used test strip and blood collection device per facility infection control policy.

Note: After a successful test, the meter turns itself off five seconds after the test strip is removed.

Flagging Test Results

It is possible to flag a test result with an asterisk (*) for a special result. When the results in memory are reviewed, this flag is a reminder of a special result.

1. Perform a test.
2. With the test result on the display, and the test strip STILL IN THE METER, press  or  once. The * appears in the lower right corner.
3. Remove and discard the used test strip and blood collection device per facility infection control policy.

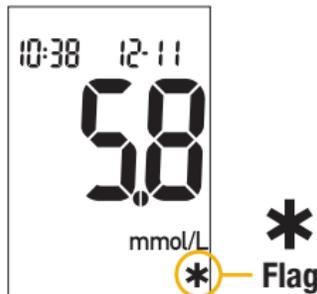Unusual Test Results

20

Follow facility policy for reporting or treating blood glucose results.

If the blood glucose result does not reflect the patient's clinical symptoms, or seems unusually high or low, run a control test. If the control test confirms that the system is working properly, repeat the blood glucose test. If the repeated blood glucose result still seems unusual, follow facility guidelines for further action or refer to the following table for assistance in problem solving.

Troubleshooting Checks

Action

1. *Check if the test strips were expired.*

Discard the test strips if they are past the expiry date.

2. *Check if the cap on the test strip container was always closed tightly.*

Replace the test strips if you think they were uncapped for some time.

3. *Check if the test strip was out of the container for a long time.*

Repeat the test with a new test strip.

4. *Check if the test strips were stored in a cool, dry place.*

Repeat the test with a properly stored test strip.

Troubleshooting Checks

Action

5. *Check if the testing steps were followed.*

Read Chapter 2, “Patient Testing for Healthcare Professionals” and test again. If you still have problems, call Roche.

6. *Check if the code number on the meter display matches the code number on the test strip container.*

If they do not match, insert the correct code chip into the meter and test again.

7. *If the problem does not resolve...*

Repeat the test with a new test strip and run a control test. If the problem does not resolve, call Roche.

Recommending Alternate Site Testing to Patients

22

Decisions about whether to recommend alternate site testing (AST) should take into account the motivation and knowledge level of the patient and his or her ability to understand the considerations relative to diabetes and AST. If you are considering recommending or using AST for your patients, you need to understand that there is a potential for a significant difference between fingertip and alternate site blood glucose test results. The difference in capillary bed concentration and blood perfusion throughout the body can lead to sample site-to-site differences in glucose results. These physiological effects vary between individuals and can vary within a single individual based upon his or her behavior and relative physical condition. Our studies involving AST of adults with diabetes show that most persons will find their glucose level changes more quickly in the fingers' blood than the alternate sites' blood.* This is especially important when glucose levels are falling or rising rapidly. If you or your patient is used to making treatment decisions based upon fingertip readings, you should consider the delay, or lag-time, affecting the reading obtained from an alternate site.

* Data on file

Chapter 3: Control Testing for Healthcare Professionals

Introduction

Testing control solutions with known glucose levels establishes that the operator and the system are performing acceptably. Control results must be within the defined acceptable ranges before valid patient testing is allowed.

Testing Intervals

Follow your facility's policy for control testing intervals.

Control testing should be performed:

- The first time before using the meter for patient testing
- At intervals established by the facility
- When a new box of test strips is opened
- If the test strip container is left open
- If test strips were incorrectly stored
- If there is a question about a patient's glucose result
- To check the performance of the system
- If the meter was dropped

Your facility may require that control testing be successfully performed **after** any of the following occur and **before** patient testing resumes:

- Previous control test results were out of the acceptable range.
- Control tests were not performed at the proper interval.

24

Control results must be within the designated range on the test strip container label, or as defined by your facility, before being considered acceptable. Patients can be tested after controls have been acceptably performed at the proper testing interval.

About the Control Solutions

- Only use Accu-Chek Performa Control Solutions.
- Write the date the bottle was opened on the bottle label. The control solution is stable for three months from that date or until the expiry date on the bottle label, whichever comes first.
- Do not use control solution that is past the expiry date.
- Close the bottle tightly after use.
- Store the bottle at 2 °C to 32 °C. Do not freeze.
- The meter automatically recognizes the Accu-Chek Performa Control Solution.
- The solution can stain fabric. Wash with soap and water.
- Control results are not displayed in memory.

Performing a Control Test

The meter, a test strip, and control solution are needed.

1. Insert the test strip into the meter in the direction of the arrows. The meter turns on.
2. Make sure the code number on the display matches the code number on the test strip container. If the code number is overlooked, remove the test strip and reinsert it into the meter.

Note: Cap color may vary – refer to the control package insert.

26

OR


(example)

3. Select the control solution to be tested. The level is entered later.

4. Put the meter on a flat surface, such as a table.

5. Remove the control bottle cap. Wipe the tip of the bottle with a lint-free wipe.

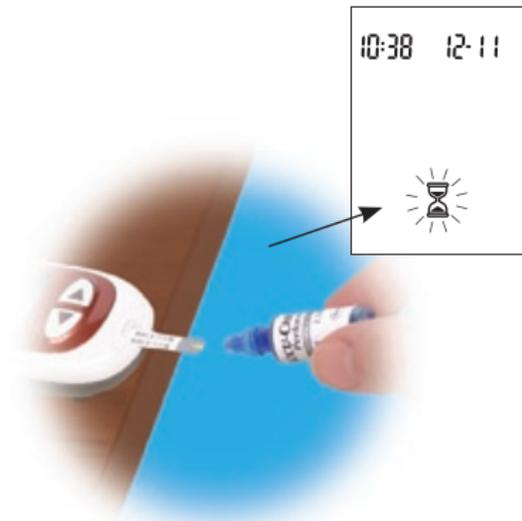
(example)

6. Squeeze the bottle until a tiny drop forms at the tip. Touch the drop to the **front edge** of the yellow window of the test strip. When  flashes, sufficient control solution is in the test strip. Wipe the tip of the bottle with a lint-free wipe, then cap the bottle tightly.
7. The result appears on the display, along with the control bottle symbol and flashing "L." Do not remove the test strip yet. Press  once to mark the result as a Level 1. Press  a second time to mark the result as a Level 2.

(example)

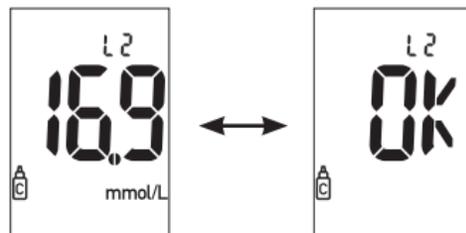
(example)

8. Press and release  to set the control level in the meter.
9. “OK” and the control result alternate on the display if the result is in range. The range is printed on the test strip

container label. “Err” and the control result alternate on the display if the result is out of range. Remove and discard the used test strip per facility policy.

Control Test Results

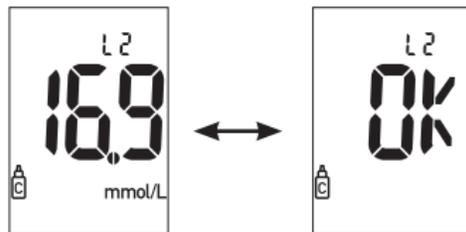(example)

The label on the test strip container shows the acceptable ranges for the Level 1 and Level 2 control solutions. Control results must be within the designated range on

the test strip container label, or as defined by the facility, before being considered acceptable.

If the control result is not acceptable, review this list for a possible cause.

1. Were the test strips or control solutions expired?
2. Was the tip of the control solution bottle wiped before and after use?
3. Were the test strip container and control solution bottle caps always closed tightly?
4. Was the test strip out of the container for a long time?
5. Were the test strips and control solutions stored properly?
6. Were the testing steps followed?
7. Was the correct control solution level selected when the test was performed?
8. Did the code number on the meter display match the code number on the test strip container?

Chapter 4: Linearity Testing for Healthcare Professionals

Introduction

Testing linearity solutions with varying glucose levels establishes the range of values from lowest to highest values for which the meter is proven capable of giving accurate results.

Follow the guidelines of your facility's regulatory agency.

Test linearity solutions the same way a control solution is tested.

Testing Intervals

Follow your facility's policy for linearity testing intervals.

It is recommended that the linearity of the system be verified the first time before using the meter for patient testing.

Linearity testing can be performed periodically to check the performance of the system

Performing Linearity Tests

The meter, test strips, and linearity solution are needed.

1. Insert the test strip into the meter in the direction of the arrows. The meter turns on.
2. Make sure the code number on the display matches the code number on the test strip container. If the code number is overlooked, remove the test strip and reinsert it into the meter.
3. Select the linearity solution to be tested.
4. Put the meter on a flat surface, such as a table.
5. Remove the linearity solution bottle cap. Wipe the tip of the bottle with a lint-free wipe.
6. Squeeze the bottle until a tiny drop forms at the tip. Touch the drop to the **front edge** of the yellow window of the test strip. When  flashes, sufficient linearity solution is in the test strip. Wipe the tip of the bottle with a lint-free wipe then cap the bottle tightly.
7. The result appears on the display, along with the control bottle symbol and flashing “L.” Do not mark the level. Record the result.
8. Turn the meter off. Remove and discard the used test strip per facility policy.
9. Repeat from step 1.

Linearity Test Results

The target values are printed in the linearity kit package insert.

Chapter 5: Meter Memory, Setup, and Downloading

Memory

Storing Control and Blood Glucose Results

The meter automatically stores up to 500 blood glucose test results with the time and date of the test. Results can be reviewed at any time. They are stored from the newest to the oldest, so set the time and date correctly in the meter.

Notes:

- The memory is not lost when you replace the battery; however, confirm that the time and date are still correct.
- Once 20 control or 500 blood glucose results are in memory, adding a new result causes the oldest one to be deleted.
- Hold  or  down to scroll through the results faster.
- The control results do not appear in memory.
- Up to 20 control results are stored in memory, but cannot be reviewed on the meter. The stored results must first be downloaded to a compatible software application. For product availability, please contact Roche.

Do not change therapy based on one individual result in memory.


Viewing Results

Press and release either ◀ or ▶ to enter memory. The most recent result appears.

Press and release ◀ to view previous results in order.

OR, press and release ▶ to look at 7, 14, or 30 day averages, in that order.

36

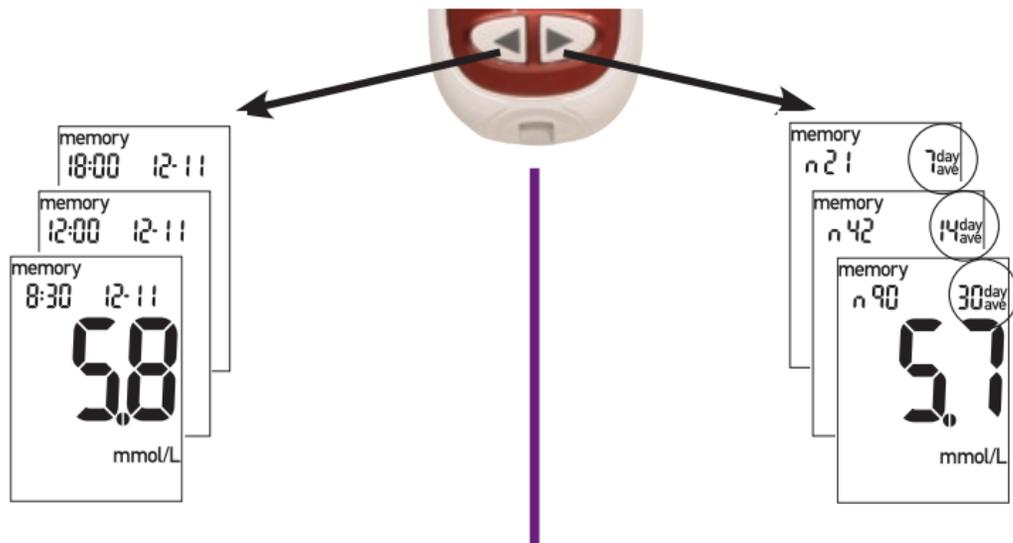Meter Setup

Using the Set-up Mode

The following features can be adjusted as required.

Time and Date – set the time and date.

Beeper – select “On” or “OFF.”

Alarm clock – select “On” or “OFF.” If “On” is chosen, select 1 to 4 alarms per day.

Hypo alarm – select “On” or “OFF.” If “On” is chosen, select the blood glucose level for the alarm.

The  has three functions for the set-up mode.

- To enter the set-up mode, turn the meter on then press and **hold**  until “set-up” flashes on the display. This takes about four seconds.
- Press and release  to set the selected option.
- To exit the set-up mode at any time, press and **hold**  for about four seconds until the flashing test strip symbol appears on the display.

You are here

TIME / DATE

(hrs, min, day,
month, year)

BEEPER

(on/off)

ALARM CLOCK


(A-1, A-2, A-3, A-4)

HYPO ALARM

(off, on, level)


Setting the Time and Date

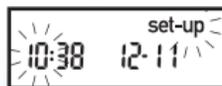38

1. Press and release  to turn the meter on. The flashing test strip symbol appears.
2. Press and **hold**  for about four seconds to enter the set-up mode. "Set-up" and the hour flash on the display.


Note: Hold  or  down to scroll faster.

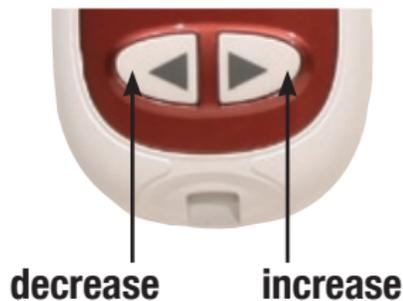
3. Press and release ◀ or ▶ to decrease or increase the hour.

4. Press and release Ⓞ to set the hour. The minutes flash.

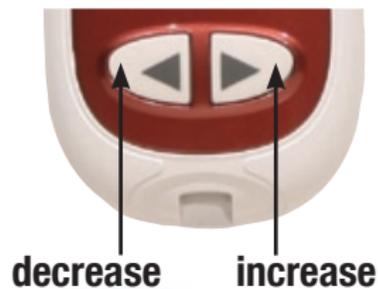5. Press and release ◀ or ▶ to adjust the minutes. Press and release Ⓞ to set the minutes.

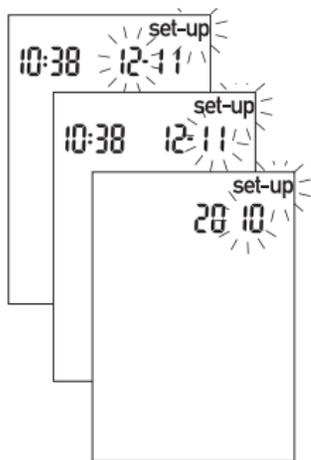
40

6. Repeat to set the day, month, and year.

7. To set up more options, press and release . To exit, press and **hold**  until the flashing test strip symbol appears on the display.

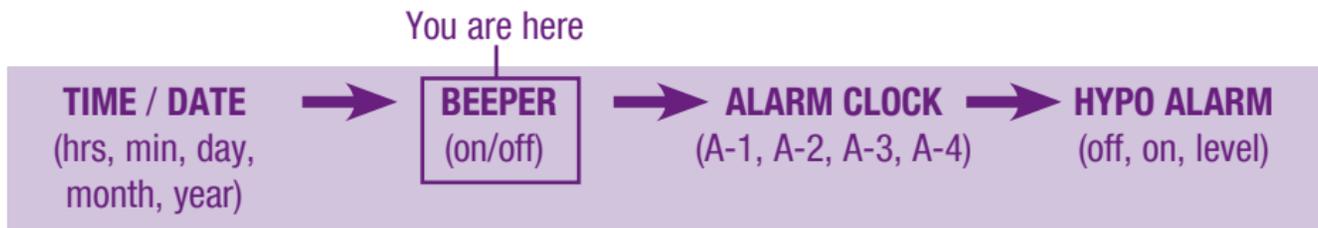Setting the Beeper On/Off


The meter beeper is preset to “On.” Setting the beeper to “OFF” does not affect test results.

The beeper prompts the operator:


- To apply blood or control solution to the test strip
- When sufficient blood or control solution is drawn into the test strip
- When the test is complete
- When a button is pressed
- If an error occurred while testing (even if the beeper is off, it beeps for an error)

1. Press and release  to turn the meter on. The flashing test strip symbol appears.

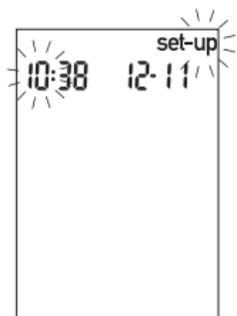2. Press and **hold**  for about four seconds to enter the set-up mode. “Set-up” and the hour flash on the display.

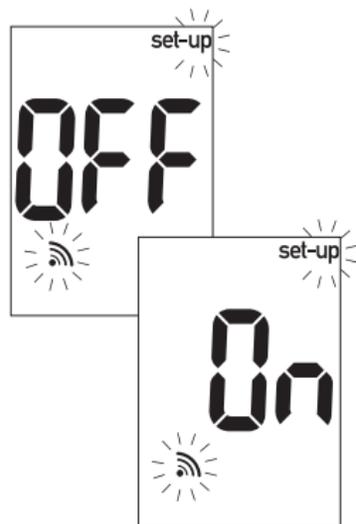
3. Press and release  repeatedly until the flashing beeper symbol and “On” or “OFF” appear on the display.

4. Press and release ◀ or ▶ to switch between “On” and “OFF.”

5. To set up more options, press and release ⓪. To exit, press and **hold** ⓪ until the flashing test strip symbol appears on the display.

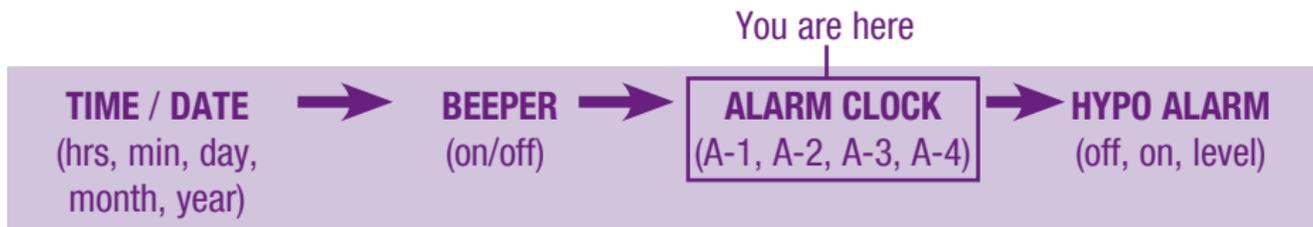Setting the Alarm Clock Function

One to four alarms can be set per day. The meter beeps every two minutes – up to three times. Turn the alarm off by inserting a test strip or pressing any button. The beeper must be set to “On” for the alarm to sound.

44

The meter is preset with the alarm clock function set to “OFF.” It must be turned “On” to use this feature.

If A-1, A-2, A-3, and A-4 are turned on, the meter is preset with the following times which can be adjusted as necessary.

A-1 8:00

A-2 12:00


A-3 18:00

A-4 22:00



Notes:

- If a test was performed within 30 minutes of an alarm, the alarm does not occur.
- If the meter is on at the alarm time, the alarm does not occur.
- Exposure to cold conditions may disable alarms until the meter is turned on.

When setting the time for the alarm clock function, the bell symbol remains on the display and “set-up” flashes continuously.


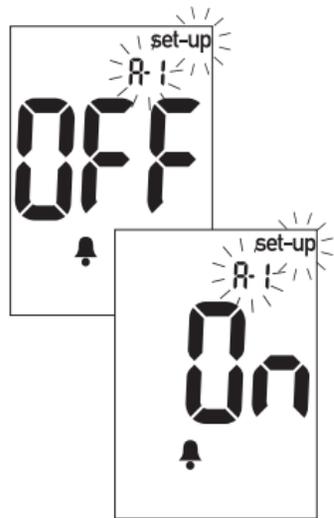
1. Press and release  to turn the meter on. The flashing test strip symbol appears.
2. Press and **hold**  for about four seconds to enter the set-up mode. “Set-up” and the hour flash on the display.
3. Press and release  repeatedly until the bell symbol, “OFF,” and the flashing “set-up” and “A-1” appear on the display.

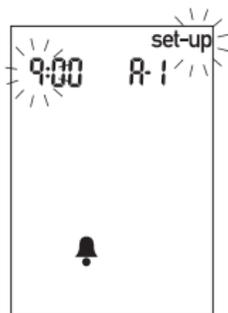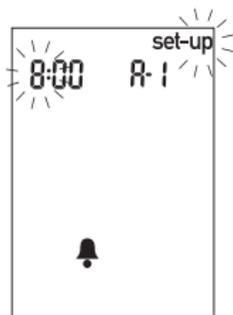- 4.** Press and release  or  to switch between “On” and “OFF.” Press and release  to set the option.

- 5.** If “On” is selected, the hour flashes. “A-1” and the bell symbol remain on the display.

- 6.** Press and release  or  to select the hour. Press and release  to set the hour.

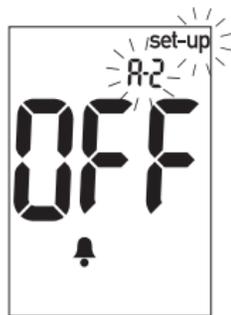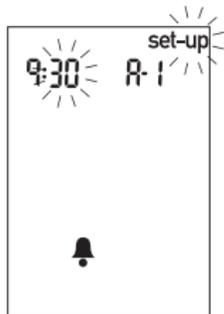7. The minutes flash. Press and release ◀ or ▶ to select 00, 15, 30, or 45. These are the only choices.

8. Press and release Ⓞ to set the minutes.

9. The next alarm “A-2” and “set-up” flash on the display with “OFF” and the bell symbol. A second

alarm can be set, or press and **hold** Ⓞ until the flashing test strip symbol appears on the display to exit the set-up mode.


Setting the Hypoglycemic (Hypo) Alarm Function

The meter hypoglycemic alarm is preset to “OFF.” The alarm level can be set from 3.3 to 4.4 mmol/L.

48

Healthcare professionals: The hypoglycemic level may vary from person to person. It is recommended to turn the hypo alarm function **off** when using the meter in a professional facility.



Consumers: Before you set the hypoglycemic alarm, talk to your healthcare professional to help you decide what blood glucose level is your hypo level.

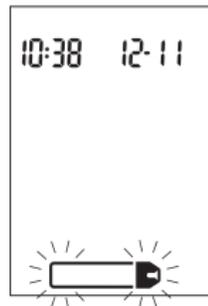1. Press and release  to turn the meter on. The flashing test strip symbol appears.




This function is no substitute for hypoglycemia training by your healthcare professional.

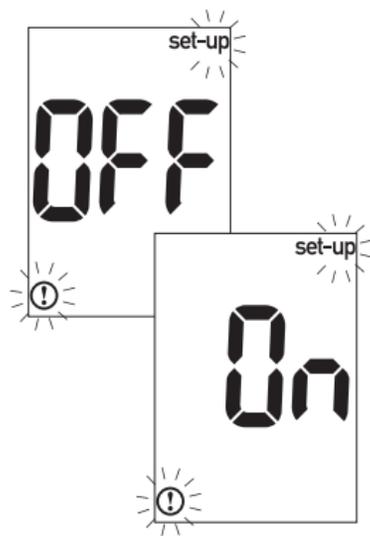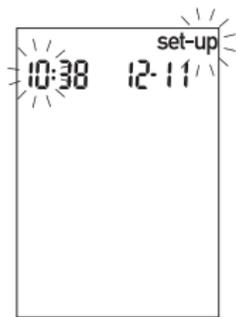2. Press and **hold**  for about four seconds to enter the set-up mode. "Set-up" and the hour flash on the display.

3. Press and release  repeatedly until "OFF" and the flashing "set-up" and  appear on the display.

4. Press and release  or  to switch between "On" and "OFF." Press and release  to set the option.

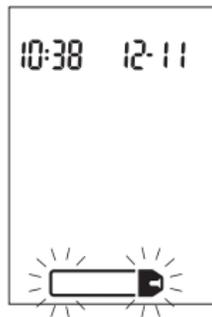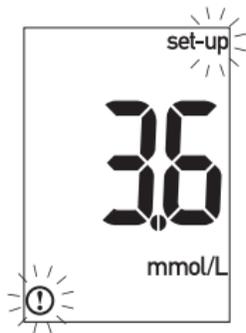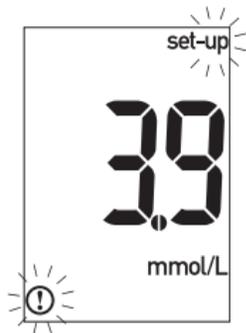5. If “On” is selected, “set-up” and  flash. The display shows 3.9 mmol/L (first use) or the current level.

6. Press and release  or  to select the level. Press and release  to set the level.


7. To exit, press and **hold**  until the flashing test strip symbol appears on the display.

Transferring Data Directly to a Computer or PDA Using Specialized Software and Infrared Cable


Downloading data is not recommended for healthcare facilities since there is no means to identify patient results.

1. Install the software according to the instructions. To transfer the results to a computer, connect the infrared meter cable according to the instructions.
2. Run the software program and follow the instructions about how to download information. Make sure the software is ready to accept data from the meter.
3. With the meter off, press and **hold** both  and  until two arrows on the display alternately flash.

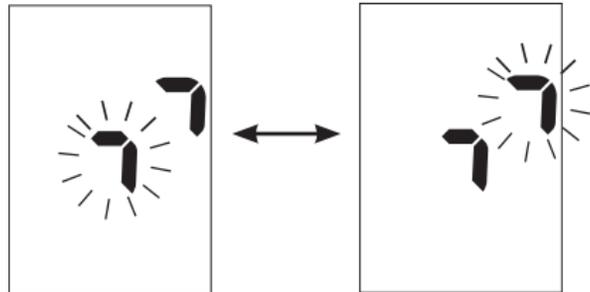4. Locate the infrared (IR) window on the top of the meter.
5. Locate the IR window on either the infrared cable (computer) or PDA.
6. Put the meter on a flat surface. Point the two IR windows toward each other. They should be 3 to 10 cm apart.
7. Do not move the infrared cable (computer), PDA, or meter during the transfer.
8. Follow the prompts on the software.
9. The software program may shut off the meter automatically when the data transfer is complete.

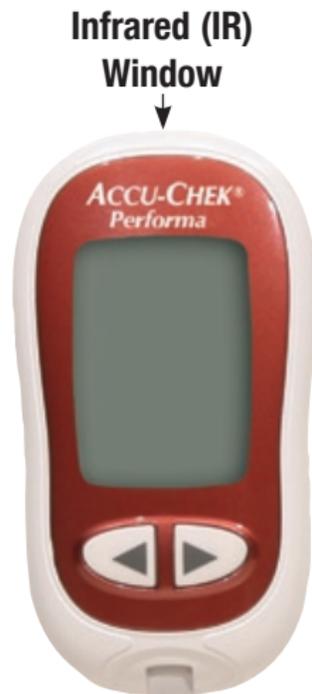Notes:

- If the data did not transfer successfully, try again. Contact Roche if problems persist.
- Be sure to set the time and date in the meter for traceability of results.

Chapter 6: Maintenance and Troubleshooting

Changing the Battery

1. Open the battery door on the back of the meter by pushing the tab in the direction of the arrow and pulling the door up. Remove the old battery.
2. Insert the new battery with the + side up.
3. Put the battery door back in place and snap it closed.
4. Turn the meter on. Confirm that the time and date are correct.

Notes:

- The meter uses one 3-volt lithium battery, type 2032.
- All control and blood glucose test results remain saved in memory.

Cleaning the Meter

Follow these instructions to clean or disinfect the meter.

- Make sure the meter is off
- Gently wipe the meter's surface with a soft cloth slightly dampened (wring out any excess) with one of these cleaning solutions:
 - Super Sani-Cloth
 - 70 % isopropyl alcohol
 - Mild dishwashing liquid mixed with water
 - 10 % household bleach solution (1 part bleach plus 9 parts water) made the same day

DO NOT

- **Get any moisture in the code chip and test strip slots**
- **Spray any cleaning solution directly onto the meter**
- **Put the meter in any type of liquid**
- **Pour liquid into the meter**

Maintenance and Troubleshooting

The meter needs little or no maintenance with normal use. It automatically tests its own systems every time it is turned on and prompts the operator for any errors. See Chapter 6, “Screen Messages and Troubleshooting.”

If the meter is dropped or the accuracy is in question, call Roche.

To check the display, turn the meter off, then press and **hold**  to see the complete display. All the segments should be clear and look exactly like the picture below. If not, call Roche.


Screen Messages and Troubleshooting


Never make treatment decisions based on an error message. If you have any concerns, call Roche.

56

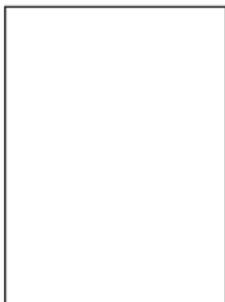
The meter will not turn on or the display is blank.

- Battery is dead –
Insert new battery
- Display is damaged –
Call Roche
- Meter is defective –
Call Roche
- Extreme temperatures –
Move the meter to a more temperate area

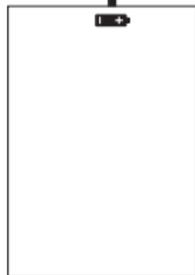Battery power is low.
Change the battery soon.

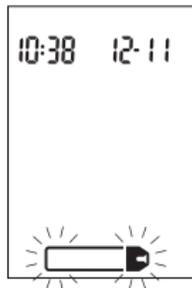The meter is ready for you to
insert a test strip.

set-up


The meter is in set-up
mode, waiting for you
to change or confirm
settings.

The meter is ready for a drop
of blood or control solution.


Blood glucose may be higher than the measuring range of the system.

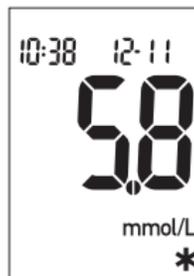This test result was flagged. See Chapter 2 or 8, “Flagging Test Results” for more details.


Blood glucose may be lower than the measuring range of the system.

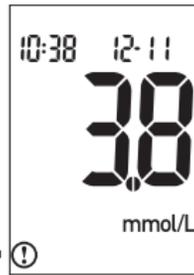Blood glucose is below the defined hypo (low blood glucose) level.

The meter is not coded or the code chip is not inserted. Turn off the meter and recode it.

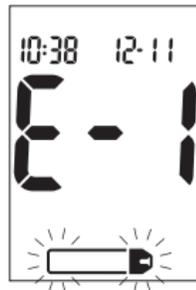The test strip is damaged. Verify that the code number on the test strip container matches the code number printed on the code chip currently inserted in the meter. Remove the test strip and reinsert it, or replace it if damaged. If the message reappears, call Roche.

The test strips will expire at the end of the current month. Before the end of the month, insert a new code chip from a new box of test strips and ensure the code chip number matches the code number on the test strip container. Make sure the time and date in the meter are correct.

The code chip is incorrect. Turn off the meter and insert a new code chip. If this does not fix the problem, call Roche.

The blood glucose may be extremely high or a meter or test strip error has occurred. If this confirms how you feel (symptoms), contact your healthcare professional immediately; healthcare professionals should follow their facility policy for treating questionable blood glucose results. If it does not confirm symptoms, repeat the test and refer to Chapter 2 or 8, “Unusual Test Results.” If this does not confirm symptoms, run a control test with the control solution and a new test strip. If the control result is within the acceptable range, review the proper testing procedure and repeat the blood glucose test with a new test strip. If the E-3 code still appears for the blood glucose test, the blood glucose result may be extremely high and above the system’s reading range. **Contact your healthcare professional immediately; healthcare professionals should follow facility policy.** If the control result is not within the acceptable range, see Chapter 3 or 9, “Control Test Results.”

Not enough blood or control solution was drawn into the test strip for measurement or was applied after the test had started. Discard the test strip and repeat the test.

Blood or control solution was applied to the test strip before the flashing drop symbol appeared on the display. Discard the test strip and repeat the test.

The code chip is from an expired lot of test strips. Ensure the code chip number matches the code number on the test strip container. Make sure the time and date in the meter are correct.

An electronic error occurred or, in rare cases, a used test strip was removed and reinserted. Turn the meter off and on, or take the battery out for 20 seconds and reinsert it. Perform a blood glucose or control test. If the problem persists, call Roche.

The temperature is above or below the proper range for the meter. Refer to the test strip package insert for system operating conditions. Move to an area with the appropriate conditions, wait five minutes and repeat the test. Do not artificially heat or cool the meter.

The battery is almost out of power. Change the battery now.

The time and date settings may be incorrect. Make sure the time and date are correct and adjust, if necessary.

Note: If you see any other error screen, please call Roche.

Chapter 7: Technical Information

Product Limitations

Please read the literature packaged with the test strips and control solutions to find the latest information on product specifications and limitations.

Specifications

Blood volume	Refer to the test strip package insert
Sample type	Refer to the test strip package insert
Measuring time	Refer to the test strip package insert
Measuring range	Refer to the test strip package insert
Test strip storage conditions	Refer to the test strip package insert
Meter storage conditions	Temperature: -25 °C to 70 °C
System operating conditions	Refer to the test strip package insert
Memory capacity	20 control and 500 blood glucose results with time and date
Automatic power off	2 minutes
Power supply	One 3-volt lithium battery (type 2032)
Display	LCD

Dimensions	93 x 52 x 22 mm (LWH)
Weight	Approx. 62 g (with battery)
Construction	Hand-held
Protection class	III
Meter type	The Accu-Chek Performa Meter is suitable for continuous operation
Control solution storage conditions	2 °C to 32 °C

Electromagnetic Compatibility

This meter meets the electromagnetic immunity requirements as per EN ISO 15197 Annex A. The chosen basis for electrostatic discharge immunity testing was basic standard IEC 61000-4-2. In addition, it meets the electromagnetic emissions requirements as per EN 61326. Its electromagnetic emission is thus low. Interference from other electrically driven equipment is not to be anticipated.

Performance Analysis

The performance data for the Accu-Chek Performa System (Accu-Chek Performa Meter with Accu-Chek Performa Test Strips) were obtained using capillary blood from diabetic patients (method comparison, accuracy), venous blood (repeatability) and control solution (reproducibility). The system is calibrated with venous blood containing various levels of

glucose. The reference values are obtained using the hexokinase method. For method comparison, the results were compared with results obtained using the hexokinase method with deproteinization (automatic analyzer). The hexokinase method is traceable to an NIST standard.

The Accu-Chek Performa System meets the EN ISO 15197 requirements.

Measuring Principle

Refer to the test strip package insert for more information.

Product Safety Information

Strong electromagnetic fields may interfere with the proper operation of the meter. Do not use this meter close to sources of strong electromagnetic radiation.

To avoid electrostatic discharge, do not use the meter in a very dry environment, especially one in which synthetic materials are present.

Disposing of the Meter

During blood glucose measurement the meter itself may come into contact with blood. Used meters therefore carry a risk of infection. Please dispose of used meters, after removing the battery, according to the regulations applicable in your country. For information about correct disposal, please contact your local council or authority.

The meter falls outside the scope of European Directive 2002/96/EC: Directive on waste electrical and electronic equipment (WEEE).

Explanation of Symbols

You may encounter the following symbols on packaging, on the type plate, and in the instructions for the Accu-Chek Performa Meter, shown here with their meaning.

Consult the instructions for use

Caution (refer to accompanying documents). Please refer to safety-related notes in the manual accompanying this instrument.

Store at

Manufacturer

REF

Catalogue number

Listed by Underwriter's Laboratories, Inc. © in accordance with UL 61010A-1 and CAN/CSA C22.2 No. 1010-1.

For in vitro diagnostic use

This product fulfils the requirements of Directive 98/79/EC on in vitro diagnostic medical devices.

3V type 2032

Guarantee

The statutory provisions on rights in consumer goods sales in the country of purchase shall apply.

Additional Supplies

The following supplies and accessories are available from authorized Roche Diabetes Healthcare Centers, pharmacies, or medical/surgical supply dealers:

Test Strips

Accu-Chek Performa Test Strips

Control Solutions

Accu-Chek Performa Control Solutions

Linearity Solutions

Accu-Chek Linearity Test Kit

Chapter 8: Blood Glucose Testing for Consumers

Why Regular Blood Glucose Testing Is Important

Testing your blood glucose regularly can make a big difference in how you manage your diabetes every day. We have made it as simple as possible.

Patients receiving peritoneal dialysis using solutions containing icodextrin (e.g. Extraneal) should not use the Accu-Chek Performa Test Strip. It could cause an overestimation of glucose results.

Performing a Blood Glucose Test

Before you perform your first blood glucose test, set up the meter correctly. You need the meter, a test strip, a lancing device, and a lancet.

1. Prepare the lancing device.
2. Wash and dry your hands.
3. Insert the test strip into the meter in the direction of the arrows. The meter turns on.

4. Make sure the code number on the display matches the code number on the test strip container. If you miss seeing the code number, remove the test strip and reinsert it into the meter.


5. A test strip and flashing blood drop symbol appear on the display.

6. Perform a fingerstick with the lancing device.


7. Gently squeeze your finger to assist the flow of blood. This helps you get a blood drop.

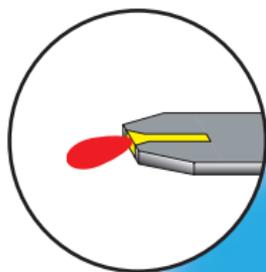10:38 12-11

8. Touch the blood drop to the **front edge** of the yellow window of the test strip. Do not put blood on top of the test strip. When you see  flash, you have enough blood in the test strip.

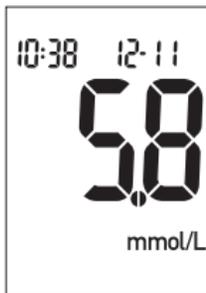9. The result appears on the display. If you want to flag the test result for a special event, leave the test strip in the meter. (See the next section.) Otherwise, remove and discard the used test strip.


Note: After a successful test, the meter turns itself off five seconds after the test strip is removed.

Flagging Test Results

If you wish, you can “flag” your test result with an asterisk (*****) to mark a special result. You might want to flag an event such as an AST result or exercise. When you review the results in memory, this “flag” can help you remember what was different about the result.


Here is how to flag a result:

1. Perform a test.
2. With the test result on the display, and the test strip STILL IN THE METER, press  or  once. The ***** appears in the lower right corner.
3. Remove and discard the used test strip.

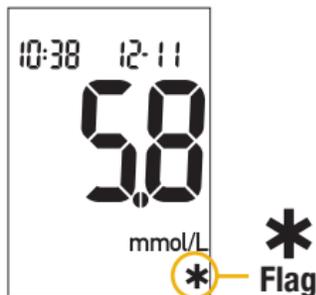Alternate Site Testing (AST)

You have the option of testing other sites on your body besides the fingertip. Blood obtained from a fingertip can be used at any time to measure blood glucose. If blood from an alternate site such as the palm, forearm, upper arm, thigh, or calf, is used, there are certain times when testing is not appropriate (see next section). The two palm testing sites are the fleshy areas under the thumb (thenar) and under the little finger (hypothenar). This is because your glucose level changes quicker in your fingertip than in the alternate sites. These differences may cause you to make the wrong therapeutic decision, producing adverse health effects. Please read the following section before you try testing from other sites.

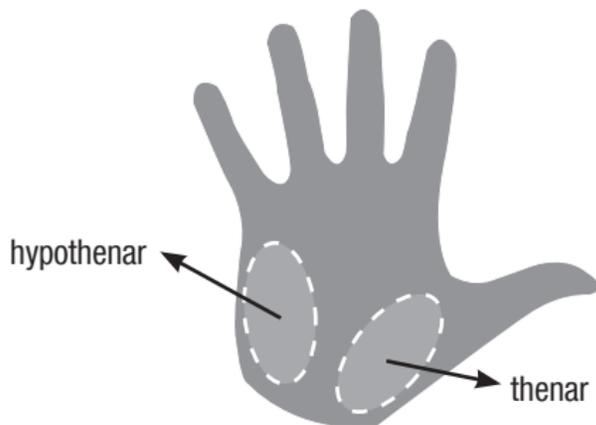IMPORTANT

- Talk to your healthcare professional about Alternate Site Testing.

Do not change your treatment because of just one result.

NEVER ignore symptoms of high or low blood glucose.

If your blood glucose result does not match how you feel, perform a fingertip test to confirm your result. If the fingertip result still does not match how you feel, call your healthcare professional.

Alternate site testing may be done:

- Immediately before a meal
- Fasting

DO NOT test from an alternate site:

- Two hours or less after eating
- After exercising
- If you are sick
- If you think your blood glucose is low
- If you often don't notice when your blood glucose is low
- During peak action time of short-acting insulin or rapid-acting insulin analogues
- Up to two hours after injecting a short-acting insulin or a rapid-acting insulin analogue

You need the meter, a test strip, a lancing device designed for AST, and a lancet.

- 1.** Prepare the lancing device.
- 2.** Insert the test strip into the meter in the direction of the arrows. The meter turns on.
- 3.** Make sure the code number on the display matches the code number on the test strip container. If you miss seeing the code number, remove the test strip and reinsert it into the meter.


4. A test strip and flashing blood drop symbol appear on the display.

5. Press the lancing device firmly against a fleshy area on the alternate site. Press the lancing device up and down in a slow pumping motion to assist the flow of blood.

6. Trigger the lancing device while keeping steady pressure on the site. Apply pressure to the site with the lancing device to assist the flow of blood.

Note: If the blood drop is too small, reapply pressure to get a sufficient blood drop.

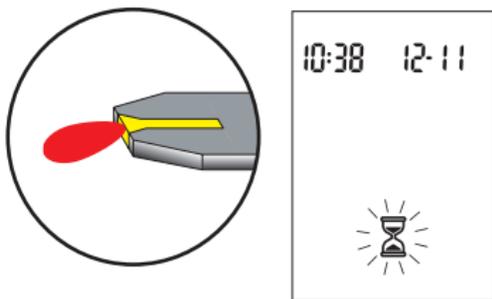7. Touch the blood drop to the **front edge** of the yellow window of the test strip. When you see  flash, you have enough blood in the test strip. Flag the result or remove and discard the used test strip.

Unusual Test Results

If your blood glucose result does not match how you feel, follow these steps:

Troubleshooting Checks

Action

1. <i>Check if the test strips were expired.</i>	Discard the test strips if they are past the expiry date.
2. <i>Check if the cap on the test strip container was always closed tightly.</i>	Replace the test strips if you think they were uncapped for some time.
3. <i>Check if the test strip was out of the container for a long time.</i>	Repeat the test with a new test strip.
4. <i>Check if the test strips were stored in a cool, dry place.</i>	Repeat the test with a properly stored test strip.

Troubleshooting Checks

Action

5. *Check if you followed the testing steps.*

Read Chapter 8, “Blood Glucose Testing for Consumers” and test again. If you still have problems, call Roche.

6. *Check if the code number on the meter display matches the code number on the test strip container.*

If they do not match, insert the correct code chip into the meter and test again.

7. *If you are still unsure of the problem...*

Repeat the test with a new test strip and run a control test. If you still have problems, call Roche.

Symptoms of High or Low Blood Glucose

Being aware of the symptoms of high or low blood glucose can help you understand your test results and decide what to do if they seem unusual. Here are the most common symptoms:

High blood glucose (hyperglycemia): fatigue, increased appetite or thirst, frequent urination, blurred vision, headache, or general aching.

Low blood glucose (hypoglycemia): sweating, trembling, blurred vision, rapid heartbeat, tingling, or numbness around mouth or fingertips.

If you are experiencing any of these symptoms, test your blood glucose. If your blood glucose result is displayed as LO or HI, contact your healthcare professional immediately.

Chapter 9: Control Testing for Consumers

Why Perform Control Tests

Performing a control test lets you know that the meter and test strips are working properly to give reliable results. You should perform a control test when:

- You open a new box of test strips
- You left the test strip container open
- You want to check the meter and test strips
- Your test strips were stored in extreme temperatures or humidity
- You dropped the meter
- Your test result does not agree with how you feel
- You want to check if you are testing correctly

About the Control Solutions

- Only use Accu-Chek Performa Control Solutions.
- The meter automatically recognizes the Accu-Chek Performa Control Solution.
- The control results are not displayed in memory.
- Write the date you opened the bottle on the bottle label. The solution is good for three months from that date or until the expiry date on the bottle label, whichever comes first.
- Do not use control solution that is past the expiry date.
- The solution can stain clothing. If you spill it, wash your clothes with soap and water.
- Close the bottle tightly after each use.
- Store the bottle at 2 °C to 32 °C. Do not freeze.

Performing a Control Test

You need the meter, a test strip, and control solution Level 1 or Level 2. The control level is printed on the bottle label.

1. Insert the test strip into the meter in the direction of the arrows. The meter turns on.
2. Make sure the code number on the display matches the code number on the test strip container. If you miss seeing the code number, remove the test strip and reinsert it into the meter.

Note: Cap color may vary – refer to the control package insert.


OR

(example)

3. Select the control solution you want to test. You will enter the level later in the test.
4. Put the meter on a flat surface, like a table.
5. Remove the control bottle cap. Wipe the tip of the bottle with a tissue.

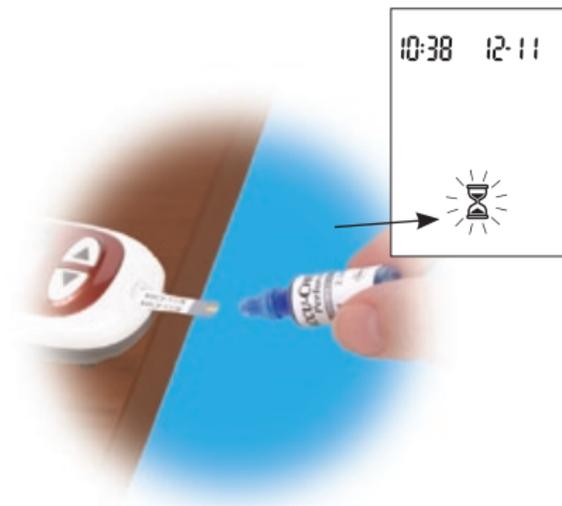
(example)

6. Squeeze the bottle until a tiny drop forms at the tip. Touch the drop to the **front edge** of the yellow window of the test strip. When you see  flash, you have enough control solution in the test strip. Wipe the tip of the bottle with a tissue, then cap the bottle tightly.
7. The result appears on the display, along with the control bottle symbol and flashing "L." Do not remove the test strip yet. Press  once to mark the result as a Level 1. If you tested the Level 2 control, press  a second time.

(example)


8. Press and release  to set the control level in the meter.
9. “OK” and the control result alternate on the display if the result is in range. The range is printed on the test strip container label. “Err” and the control

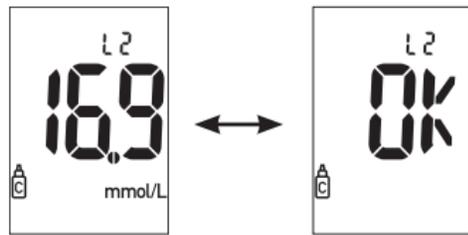(example)

result alternate on the display if the result is not in range. Remove and discard the used test strip.

Understanding Control Test Results

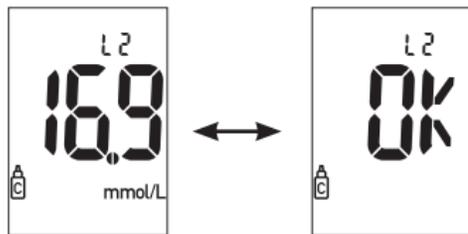(example)

The label on the test strip container shows the acceptable ranges for both the Level 1 and Level 2 control solutions. The result you get should be inside this range. Make sure you compare the result to the correct level of control.

When the control result is inside the range on the test strip container, the test strips and meter are working properly.

If the control result is not inside the acceptable range, here are some things you can do to solve the problem:

Troubleshooting Checks

Action

1. *Check if the test strips or control solutions were expired.*

Discard the test strips or control solution if either is past the expiry date. If the control solution was opened more than three months ago, discard it.

2. *Check if you wiped the tip of the control solution bottle before and after use.*

Wipe the tip of the bottle with a tissue. Repeat the control test with a new test strip and a fresh drop of control solution.

3. *Check if the caps on the test strip container and the control solution bottle were always closed tightly.*

Replace the test strips or control solutions if you think either was uncapped for some time.

Troubleshooting Checks

Action

4. Check if the test strip was out of the container for a long time.

Repeat the control test with a new test strip.

5. Check if the test strips and control solutions were stored in a cool, dry place.

Repeat the control test with a properly stored test strip or control solution.

6. Check if you followed the testing steps.

Read Chapter 9 “Control Testing for Consumers” and test again. If you still have problems, call Roche.

7. Check if you chose the correct control solution level, either 1 or 2, when you performed the test.

If you chose the wrong control solution level, you can still compare the control result to the range printed on the test strip container.

Troubleshooting Checks

Action

8. *Check if the code number on the meter display matches the code number on the test strip container.*

If they do not match, insert the correct code chip into the meter and test again.

9. *If you are still unsure of the problem...*

Repeat the control test with a new test strip. If you still have problems, call Roche.

Index

- alarm clock function,
setting, 44
- alternate site testing,
Consumers, 75
Healthcare professionals, 22
- appetite, increased, 82
- battery, changing, 53
- battery, installing, 53
- battery, type, 54
- beeper, setting, 41
- blood glucose, testing,
Consumers, 70
Healthcare professionals, 15
- blurred vision, 82
- button, on/off/set, 7, 8
- code chip, 8, 9
- computer, connecting the
meter to, 51
- control solution,
Consumers, 84
Healthcare professionals, 24
- control test, acceptable
range,
Consumers, 89
Healthcare professionals, 29
- control test, performing,
Consumers, 85
Healthcare professionals, 25
- control test, results,
Consumers, 89
Healthcare professionals, 29
- control test results,
understanding,
Consumers, 89
Healthcare professionals, 29
- control test, unacceptable
range,
Consumers, 90
Healthcare professionals, 30
- display check, 55
- error messages, 56–62
- expiry date, 13, 24, 84
- fatigue, 82
- flagging test results,
Consumers, 74
Healthcare professionals, 19
- frequent urination, 82

- guarantee, 68
- high blood glucose,
 - Consumers, 82
 - Healthcare professionals, 20
- hyperglycemia,
 - Consumers, 82
 - Healthcare professionals, 20
- hypoglycemia,
 - Consumers, 82
 - Healthcare professionals, 20
- hypoglycemic alarm
 - function, setting, 48
- linearity, testing, 31
- low blood glucose,
 - Consumers, 82
 - Healthcare professionals, 20
- maintenance, meter, 53
- memory, meter, 35
- numbness, 82
- product limitations, 63
- product specifications,
 - 63–64
- rapid heartbeat, 82
- sample handling, healthcare professionals, 15
- settings, meter, 37
- supplies, 68
- sweating, 82
- symbols, 67
- technical information, 63
- test results, unusual,
 - Consumers, 80
 - Healthcare professionals, 20
- test strips, Accu-Chek Performa, 8, 13
- thirst, increased, 82
- time and date, setting, 38
- tingling, 82
- trembling, 82
- troubleshooting, 56–62

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ACCU-CHEK and ACCU-CHEK PERFORMA are trademarks of Roche. All other product names and trademarks are the property of their respective owners.
©2009 Roche Diagnostics. All rights reserved. / Tous droits réservés. 05104491001-1009
Made in U.S.A. for distribution in Canda by:
Roche Diagnostics GmbH